


## Processing Sinks

Our processing sinks are designed to facilitate a wide variety of cleaning tasks in a safe and comfortable manner. That's safe and efficient manual cleaning of medical instruments and devices.

Engineers of Confidence.

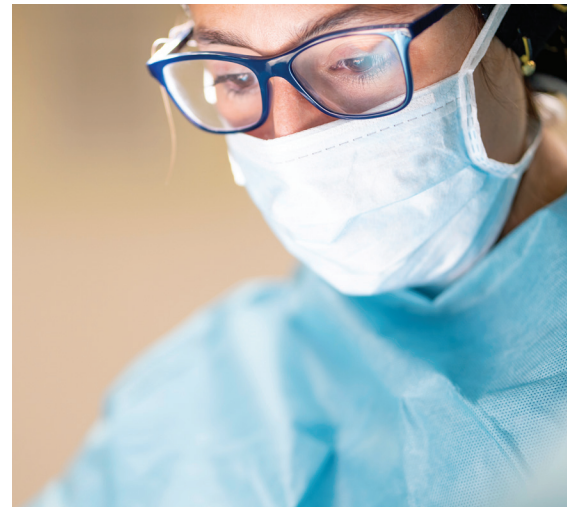

# Processing Sinks

### Versatile

We offer single, double, and triple basin sinks to fit a wide variety of applications and cleaning needs. Sinks may be configured with the desired cleaning tools and workflow direction to optimize clinician efficiency during the cleaning process. An optional pegboard, with a full portfolio of accessories, improves organization and efficiency.

### Durable

Belimed's processing sinks feature rigid, heavy-gauge stainless steel construction with seam-welded components. Our sinks are built to provide years of reliable operation without excess vibration, distortion, sagging, or other damage.

### Ergonomic

Our powered, height-adjustable models can accommodate a wide range of clinician heights. These sinks are built to the AAMI HE75 2009 and ANSE/HFS 100-1988 standards and can meet the needs of the 5th through 95th percentile of users. Memory preset controls permit one-touch countertop height adjustment and each user may quickly customize the sink surface to fit his/her stature.

### Optional Features

A variety of options and accessories are available to enhance our processing sinks' capabilities and workflow value.

#### Options:


- Drain Configuration (PVC or Copper)
- Side Splashes
- Cleaning Faucets
- DI Faucet
- Workflow Direction
- Air and Water Guns with Tips
- Ultrasonic Basin

#### Pegboard with Accessories:

- LED Overhead Light
- Glovebox Holder
- Bottle Holder
- Gallon Container Holder
- Brush Holder
- Bin Rail
- Small Bin and Large Bin
- Display Arm and Mount Bracket
- Shelf



Versatile

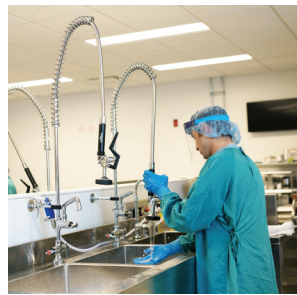

Durable



Ergonomic

Distributed by Belimed for Mac Medical

Feature / Option	Single Basin Sink	Double Basin Sink	Triple Basin Sink
Number of Basins	1	2	3
Basin Widths (in)	24 or 30 (in)	24 or 30 (in)	24 or 30 (in)
Sonic Basin option	Yes	Yes	Yes
Backsplash	12" H, Sloped Top	12" H, Sloped Top	12" H, Sloped Top
Sidesplash	10" H, Left, Right or Both	10" H, Left, Right or Both	10" H, Left, Right or Both
Faucets	Pre-rinse/Gooseneck combo or Gooseneck only	Pre-rinse/Gooseneck combo or Gooseneck only	Pre-rinse/Gooseneck combo or Gooseneck only
DI Faucet	Yes	Yes	Yes
Air and Water Guns	Yes, with Tips and Mount Bracket	Yes, with Tips and Mount Bracket	Yes, with Tips and Mount Bracket
Drain Types	Strainer or Lever	Strainer or Lever	Strainer or Lever
Fixed Height	Yes	Yes	Yes
Adjustable Height	Yes	Yes	Yes
Pegboard & Accessories	Yes	Yes	Yes
Dimensions, (in) W X D X H	72 x 30 x 37, Fixed Height (in) 72 x 30 x 34-42, Adjustable Height (in)	96 x 30 x 37, Fixed Height (in) 96 x 30 x 34-42, Adjustable Height (in)	120 x 30 x 37, Fixed Height (in) 120 x 30 x 34-42, Adjustable Height (in)

REV4.5.2021\_Q320

Engineers of Confidence.