

## Product Specification

# Single Basin Processing Sink

### Product Description

A single basin, all stainless, processing sink with options for basin size, sonic basin, faucets, guns, drain types, side splashes and fixed or adjustable height function.

### Application

Used in Sterile Processing decontamination area to perform pre-cleaning for instruments, containers and other related items.

### Standard Features

#### Construction/Design

- Type 300 stainless construction.
- 12" high integrated backsplash with sloped top and mounting holes for faucets and accessories
- Solid sink deck, with no mounted accessories, for easy cleaning
- Sloping drain board on each side
- No-drip marine edge on sink deck
- Basins are etched with water levels indicating gallons

### Dimensions

- Fixed Height Sink - 72" W x 30" D x 37" H
- Adjustable Height Sink - 72" W x 30" D x 34" to 42" H

### Basin Width Options

- 24" W x 17" D x 10" H
- 30" W x 17" D x 10" H

### Benefits

#### Rinse and Fill Faucets

Sink can be configured with up to two faucets. Pre-rinse/Gooseneck combination faucet and a gooseneck only faucet are available. All faucets mount to the backsplash. Pre-rinse/Gooseneck combo faucets also include a clamp to secure the faucet to the wall or pegboard.

#### DI Faucet

An ON/OFF knob controlled, tin-lined, DI faucet is available. Mounts to the backsplash.

#### Air and Water Guns

Air and water guns, both with coiled hoses, hooked body design, tips and wall or pegboard mounting brackets are available. Both gun connections mount to the backsplash.

#### Drain

Drain is available as standard strainer type and lever operated. PVC tailpiece included with fixed height units. Drain piping is included for adjustable height sinks and can be specified as either copper or PVC.



#### Side Splashes

10" high side splash available for right side, left side or both sides.

#### Ultrasonic Basin

The basin can be configured as an ultrasonic cleaner, for either the 24" or 30" wide basins. The ultrasonic power generator is located under the right-hand drain board. A simple knob is used for setting sonic on time. This option requires the use of lever drains.

#### Adjustable Height Control Panel

When configured as an adjustable height unit a small control panel is mounted, on the front, under the right-hand drain board. The panel controls up and down movement and allows for up to 4 user settable heights.

#### Pegboard & Accessories

A full length, stainless steel, pegboard, with ¼" holes on a 1" pattern is available. A full portfolio of accessories are available and can be added at any time.

#### Utilities / Installation Requirements

Reference cut sheet for utilities and installation requirements.



For more information, visit [belimed.com](http://belimed.com)

# Single Basin Processing Sink

### Options

#### Type

- Fixed Height
- Adjustable Height

#### Basin Widths

- 24" wide basins
- 30" wide basins

#### Basin Type

- Ultrasonic Basin

#### Drain Type

- Standard Strainer Drain
- Lever Operated

#### Drain Material (Adjustable Height)

- PVC
- Copper

#### Side Splash


- Left Side
- Right Side
- Both Sides

#### Faucet Type

- Pre-rinse/Gooseneck
- Gooseneck Only

#### DI Faucet, Air and Water Guns

- Treated (DI) faucet
- Air Gun with tips
- Water Gun with tips

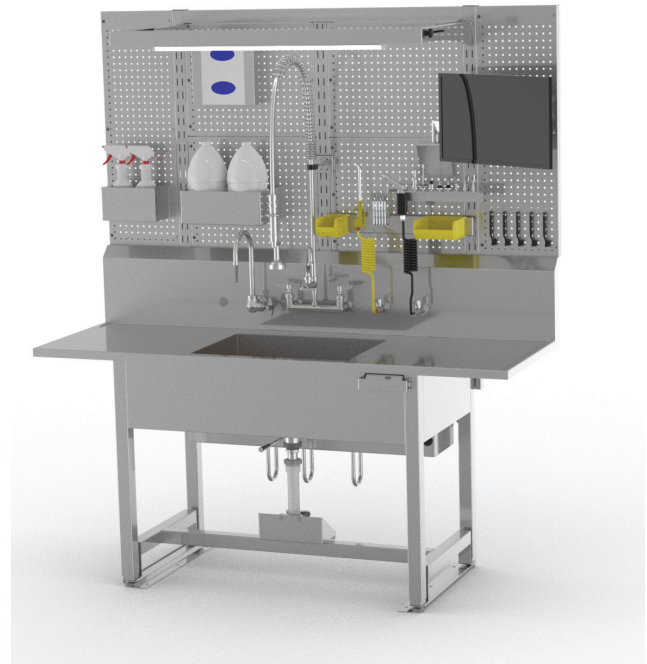

### Pegboard & Accessories Options

#### Pegboard

- 72" W x 36" H Pegboard

#### Accessories

- LED Overhead Light
- Glovebox Holder
- Bottle Holder
- Gallon Container Holder
- Brush Holder
- Bin Rail
- Small Bin
- Large Bin
- Display Arm and Mount Bracket
- Shelf

Distributed by Belimed for Mac Medical

REV4.5.2021\_Q0

Engineers of Confidence.

**Belimed**  
Infection Control

Belimed, Inc.  
1.800.451.4118 | 1.843.216.7424  
[belimed.com](http://belimed.com)