



WA 290 Automation System

**Automated Loading and Unloading for the
WD 290 IQ Washer / Disinfector**

Engineers of Confidence.

Belimed
Infection Control


WA 290 Automation System

Decontamination Side

Designed to make the decontamination side more effective and efficient, for the way your team works, with two types of loading options.

Shuttle Loading system (WA 290 CL)

- Load racks onto the queuing conveyor system.
- Shuttle picks up a rack and delivers to the next available washer, loads washer and starts the cycle.
- Shuttle can be "suspended" to manually load an emergency rack into an open washer.
- Various load queue capacities to meet your space and throughput needs.
- Provides a single point loading, from a cart, onto the system.

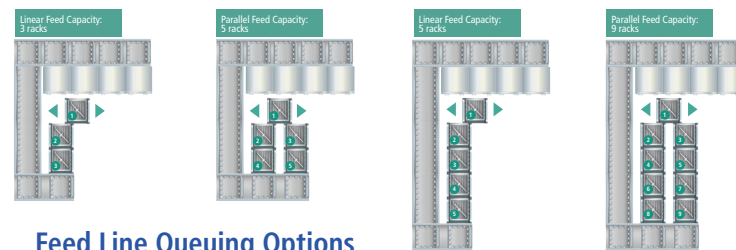


Feed Line Queuing Options



Increased Efficiency & Performance = Reduction of Workload and Costs

Dedicated Conveyor Loading (WA 290 SI)

- Load racks onto the input queuing conveyor in front of a specific washer.
- Conveyor will auto load the rack into the washer, and start the cycle when the washer becomes available.
- Single rack load conveyor capacity.
- Provides a more familiar method of direct loading from a cart.

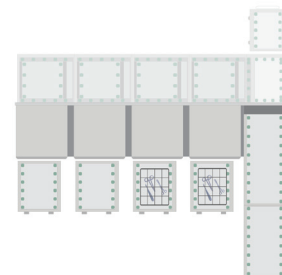

Clean Side

Designed to make the clean side more effective and efficient, for the way your team works, with two types of unloading options

Cross Conveyor Unloading System (WA 290 CL)

- Washer automatically unloads rack onto the cross-conveyor system.
- Rack travels to unload/load station in front of the return window.
- Full racks can be removed onto a cart at this location.
- Empty racks can be loaded to return to the Decontam side.
- Racks can also be left on the conveyor and clean instrument removed from the rack.
- Unload conveyor capacity of 1 rack.
- Provides a single point unloading and rack return.



Dedicated Conveyor Unloading (WA 290 SO)

- Washer automatically unloads racks onto the output conveyor.
- Dock cart and then automatically unload from the conveyor to the cart.
- Empty racks are loaded onto the return line to go back to the decontam side.
- Various return line rack capacities to fit your throughput and space needs.
- Provides a more familiar method of direct unloading to a cart.



Contact

Belimed Inc. is a leading supplier of medical and surgical instrument sterilization, disinfection, and cleaning products and services.

With 12 Belimed companies throughout Europe, North America and China, Belimed is represented worldwide by a strong network in over 80 countries.

USA Headquarters

Belimed, Inc.
8351 Palmetto Commerce Parkway
Suite 101
Ladson, SC 29456
1.800.451.4118
info.us@belimed.com

Global Headquarters

Belimed AG
Grienbachstrasse 11
6300 Zug
Switzerland
+41 41 449 78 88
info@belimed.com

REV3.8.2023_Q0



Engineers of Confidence.

Belimed
Infection Control

For more information,
visit www.belimed.com

