




Belimed Protect® Pretreatment Foamer

The Belimed Protect® Pretreatment Foamer is designed to easily and ergonomically apply Belimed Protect Pretreatment to surgical instruments. The foamer includes a unique cart for mobility and two wands, one for general trays and one for cannulated instruments.

Engineers of Confidence.

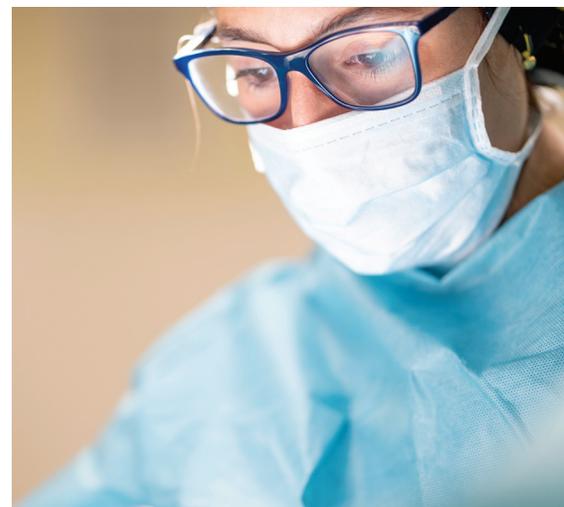

Belimed Protect® Pretreatment Foamer

Standard Features

- The Pretreatment Foamer allows for ergonomic application for staff, resulting in better compliance for pretreating instruments
- Two Application Wands – one for general instruments and one specially designed for cannulated instruments
- Triple the coverage of manual spraying in the same application time
- The foam created from the Pretreatment Foamer minimizes aerosolization
- The cart design holds 1 gallon or 2.5 gallon container sizes, which is more cost effective and environmentally friendly than traditional manual spray bottles



Specifications

- Part Number: 1064361 (Unit), 1070822 (Cart)
- Working Voltage: 110-120VAC
- Power Consumption: 30W
- Frequency: 60Hz
- External Circuit Insulation Rating:
E495378 (UL) SVT 3*18AWG(0.824MM)
105°C 300V
VW-1 C(UL) SVT 3*18AWG(0.824MM)
105°C 300V FT1
Insulation is grade A



Easy Access Storage



Solution for Cannulated Instruments




Intertek
5025042

Conforms to UL STD 61010-1
Certified to CSA STD C22.2# 61010-1-12

Distributed by

Belimed
Infection Control

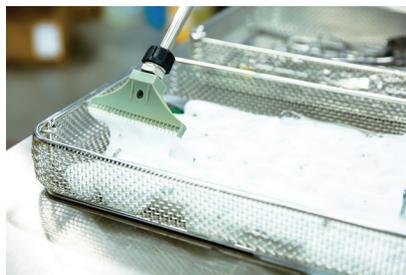

Ergonomic Application



Triple the Coverage

DEV3.16.2023_Q0

Engineers of Confidence.

Belimed
Infection Control

Belimed, Inc.
1.800.451.4118 | 1.843.216.7424
belimed.com