

## Product Specification

# Belimed Protect® Pretreatment Foamer

### Product Description


The Belimed Protect® Pretreatment Foamer is designed to easily and ergonomically apply Belimed Protect Pretreatment to surgical instruments. The foamer includes a unique cart for mobility and two wands, one for general trays and one for cannulated instruments.

### Standard Features

- The Pretreatment Foamer allows for ergonomic application for staff, resulting in better compliance for pretreating instruments
- Two Application Wands – one for general instruments and one specially designed for cannulated instruments
- Triple the coverage of manual spraying in the same application time
- The foam created from the Pretreatment Foamer minimizes aerosolization
- The cart design holds 1 gallon or 2.5 gallon container sizes, which is more cost effective and environmentally friendly than traditional manual spray bottles

### Specifications

- Part Number: 1064361 (Unit), 1070822 (Cart)
- Working Voltage: 110-120VAC
- Power Consumption: 30W
- Frequency: 60Hz
- External Circuit Insulation Rating: E495378 (UL) SVT 3\*18AW G(0.824MM) 105°C 300V, VW-1 C(UL) SVT 3\*18AWG(0.824MM) 105°C 300V FT1, Insulation is grade A

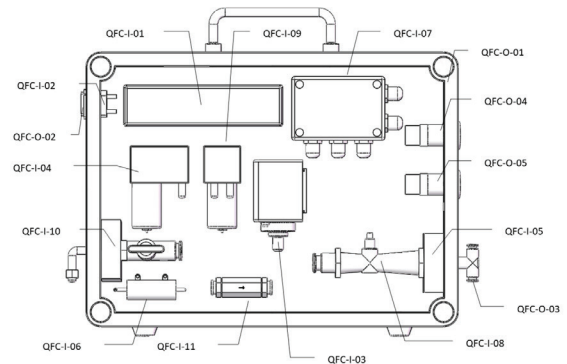

Item Number	Description
QFC-I-01	Power Supply
QFC-I-02	Fuse
QFC-I-03	Pressure Control
QFC-I-04	Air Pump
QFC-I-05	Foaming Out Port
QFC-I-06	Silencer
QFC-I-07	Control Box
QFC-I-08	Foam Generator
QFC-I-09	Water Pump
QFC-I-10	Detergent Entry
QFC-I-11	Check Valve
QFC-O-01	Cabinet
QFC-O-02	Power Outlet
QFC-O-03	Hose Connection
QFC-O-04	Power Switch
QFC-O-05	Liquid Wash



Conforms to UL STD 61010-1  
Certified to CSA STD C22.2# 61010-1-12

Distributed by  
**Belimed**  
Infection Control



For more information, visit [belimed.com](http://belimed.com)