

Product Specification

Prep and Pack Ergonomic Workstations

Designed with clinician comfort and safety in mind, our prep and pack ergonomic workstations are ideal for a variety of instrument processing tasks.

Product Description

Belimed's line of Prep & Pack workstations are designed to promote efficient, ergonomic processing and packaging of instruments and utensils. Workstations are highly configurable, available in various sizes, with different work surface materials, and with a wide array of storage, lighting, and technology support options. Frames are also available for each size, in fixed, manual adjust, and power adjustable configurations.

Application

For use in healthcare facilities for the inspection, organization, and wrapping of re-usable instruments and utensils prior to sterilization.

Standard Features

Construction/Design

Belimed Prep & Pack workstations are constructed of high-quality materials and finishes and include the following standard features:

- Type 304 stainless steel standard tops with seam welded corners and vibration-damping solid core
- Heavy-gauge steel frames with durable powder-coat paint surfaces
- Cantilevered frame design (powered & manual) for maximum knee well access
- Up to 1,200 lb. lift capacity on electric adjust units
- Slotted uprights for tool-free overhead storage reconfiguration
- Widest range of accessories

Benefits

Optimized Ergonomics

Manual height adjustment (ideal for the single shift/operator station), and electric powered height adjustment (best choice for multiple operator/shift applications) are available. Programmable controls that store user height settings are a standard option. Fixed-height is also available as an entry-level option. All workstations and accessories are designed to adjust to fit the needs of 5th through 95th percentile users, having been built to the ANSI/HFS 100-1988 standard.

Ideal for Multiple Users

Power adjust models feature smooth action with simple easy to use under-counter up/down controls. Each user may select the proper working height to fit their stature, thereby ensuring hours of fatigue-free operation. A height-adjustable workstation allows the users to adjust to an optimized, neutral position and reduce or eliminate the ergonomically risky and non-value adding body movements while working.

Massive Stability

All frames are fully welded of heavy-gauge steel and include back rails to prevent frame distortion and racking under load. Solid core stainless tops reduce vibration and sagging common to standard welded-gusset reinforced-top designs.

Easy Reconfiguration

Tool-free overhead storage system connections allow shelves and bins to be quickly and easily reconfigured, to meet changing department or user needs.

Serviceability

Motors and controls for power adjustable units drop out for quick and easy change-out with minimal downtime.



Options

Surface Material:

- Stainless Steel
- T-Mold Plastic/Cleanroom Laminate

Frames:

- Power
- Manual
- Fixed

Workstation Widths:

- 48"
- 60"
- 72"

More sizes available on request

Work Surface Depths:

- 30"
- 36"

More sizes available on request

Wide Range of Storage & Organizational Options:

Shelves, Bin racks, Drawers, Wrap rods, Rails, and Baskets are all available to keep the supplies you need close at hand.

NEW OPTION – Social Distancing Dividers:

These can be used to quickly create separate workspaces. These can form a single separation wall, an L-shaped enclosure, a U-shaped enclosure, and many more configurations. The polycarbonate panel allows users to see through and it can be easily cleaned, supporting your sanitation protocol.

Full Technology Support:

PC's, printers, and scanners slot easily into Belimed's workstation, with mounting options for flat screens, CPU's, printers, and keyboards, and comprehensive power outlet and cord management support. Auxiliary carts printer/PC carts are also available to help maximize countertop and under counter space.



For more information, visit belimed.com

Prep and Pack Ergonomic Workstations

Accessories

Overhead Storage:

- 8" x 30" shelf, 12" x 30" shelf, 15" x 30" shelf
- 4" x 30" bin rail
- Spool holder rod
- Small roll shelf
- Peel pack box dispenser (for bin rail)
- Wrap paper backframe (4 rods above work surface)

Under-counter Storage:

- Single 6" drawer, Single 12" drawer
- Dual-stacked 6" drawers
- Dual slide-out wire baskets

Power Management:

- 8 outlet power bar
- Rear wire trough
- Horizontal power bar work surface mount

PC Support:

- Articulating keyboard arm with mouse tray (fixed or adjustable height)
- Articulating flat-panel display arm (fixed or adjustable height)
- Fixed shelf mount keyboard tray
- Pull-out/Under-surface keyboard tray
- CPU undermount bracket

Overhead Lighting:

- 48" Dual-tube
- 36" Dual-tube

Task Lighting:


- Fluorescent bar light
- Fluorescent light ring magnifier
- LED light ring magnifier

Organizational Aids:


- Magnetic marker whiteboard
- Pull-out work surface
- Pegboard

Miscellaneous:

- Adjustable footrest
- Counter backplash
- Locking & non-locking casters
- Leveling casters
- Printer carts
- Wrap paper rack
- Tool balancer trolley

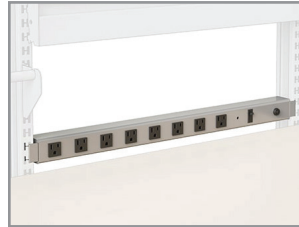


8 Outlet Power Bar

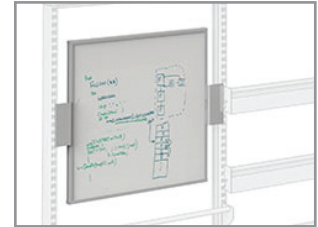


Magnetic Whiteboard



Flat Panel Display Arm

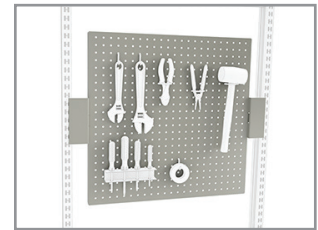

Pegboard

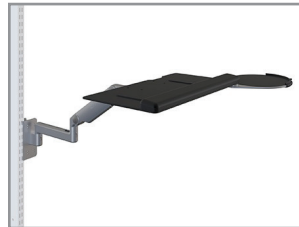

Keyboard Arm with Mouse Tray



Free Standing Footrest



Light Ring Magnifier



Tool Balancer Trolley

Distributed by Belimed for BOSTONtec

REV3.3.2023_Q0

Engineers of Confidence.

Belimed
Infection Control

Belimed, Inc.
1.800.451.4118 | 1.843.216.7424
belimed.com